



To use the plugin, download the ImageJ macro Click_cells2.ijm below and open ImageJ (Fiji).
Choose from the menu bar Plugins > Macros > Install.

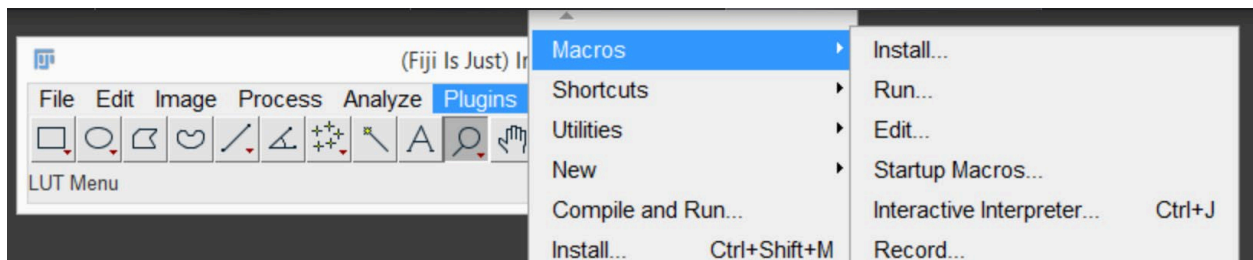

Choose Click_cells2.ijm


A new button will appear in the ImageJ toolbar

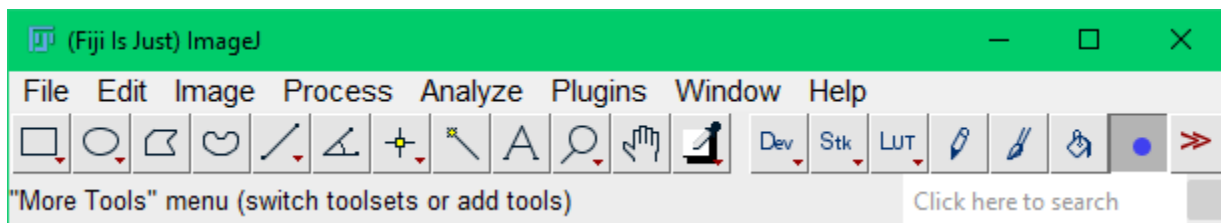

Open the measurement movie that you want to annotate.

Click on the new tool (round blue dot) to activate the annotation tool.

Annotation can be done by clicking on the cell that you want to annotate. A new row will be added to a table.

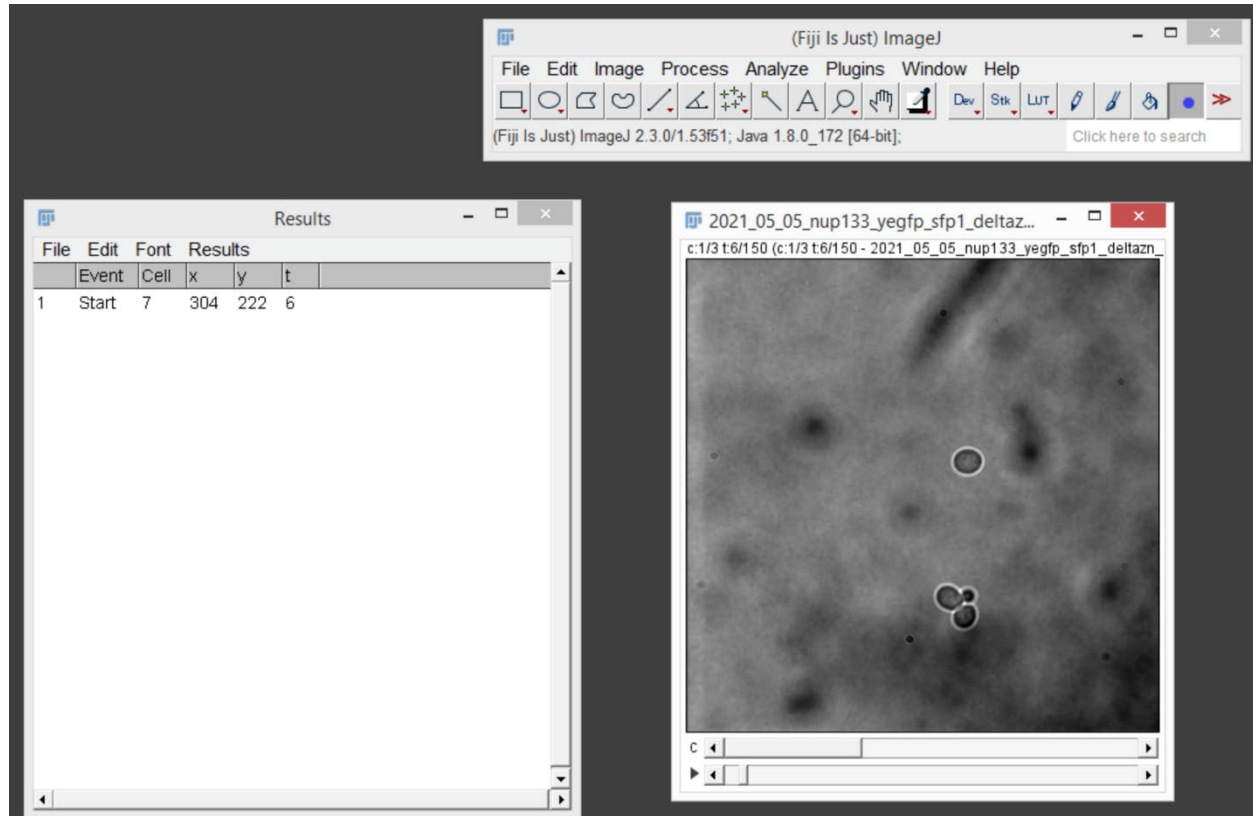

Clicking somewhere on the image will start a new cell (numbering goes automatically; it increases the cell number). You can now add other events such as budding, cytokinesis, cell death, cell wash. You can annotate by using the following keys:

- a - Cell budding
- c - Cytokinesis
- d - Cell death
- w - Cell washed

Other functionality:

- n - Increase cell number
- p - Decrease cell number

s - Settings; this can be used to set the number of channels (default is 3). This is important since it is used to calculate the correct frame number.

e - Export the results to a format that can be read by the Python analysis program.

If you made a mistake, you can always adjust the results table by deleting rows. Select the row you want to delete and then select Edit > Clear.

When you have annotated all cells, you can save the results table and/or export everything to a text file by pressing 'e'. When you have a results' table you can reopen that table later and continue annotating with this tool.